

11.6" Win10 Tablet Digital Camera, 16.0M



- 2.4G/5G WIFI, Bluetooth, USB3.0, HDMI, Ethernet Port
- S-EYE Measure Software Installed
- Type C Power Supply DC12V/2A Adapter

Product Details



11.6" LCD, Resolution 1920x1080p IPS Touch Screen Win10 Tablet





16.0 M Digital Camera All-in-One Design

Product Details



Software For Windows

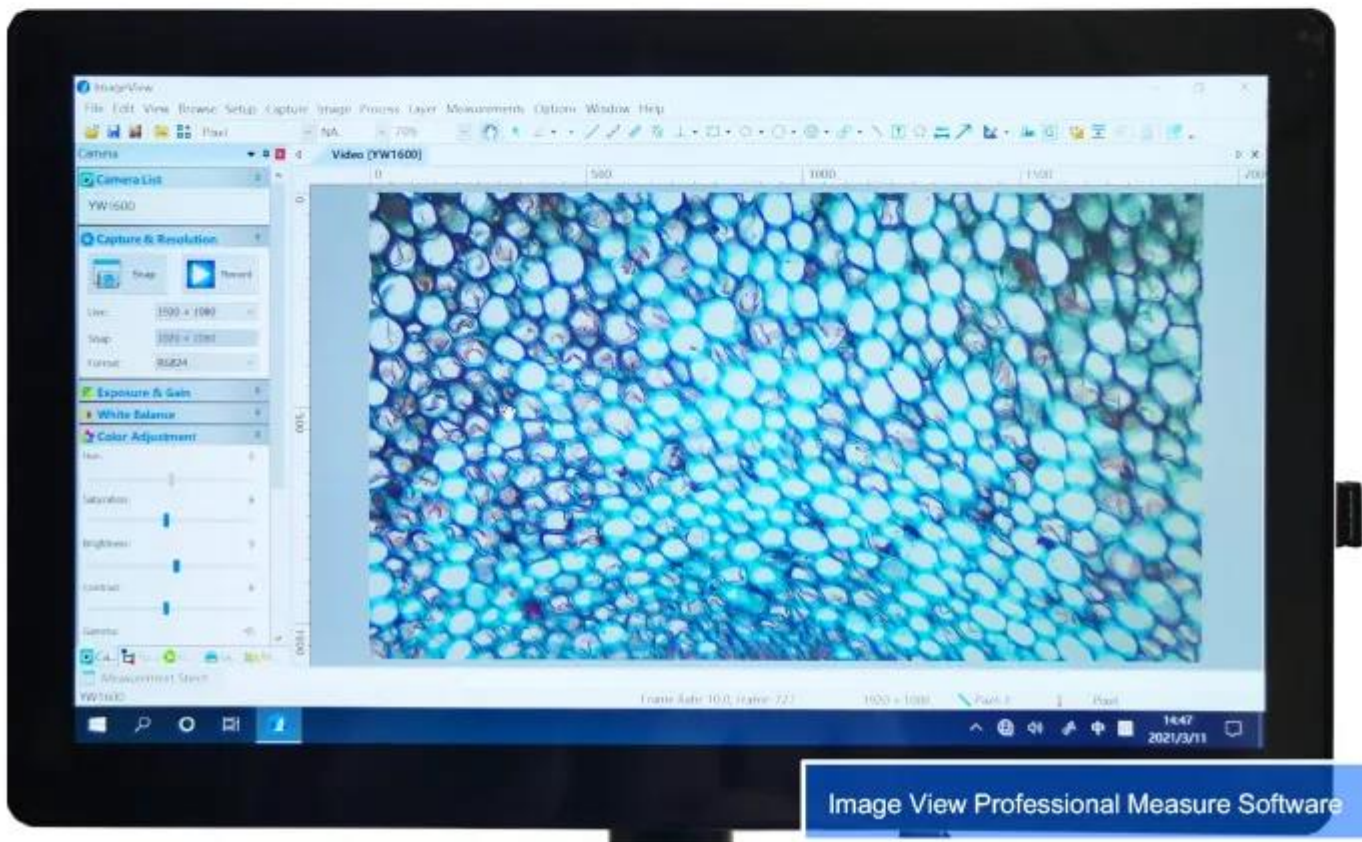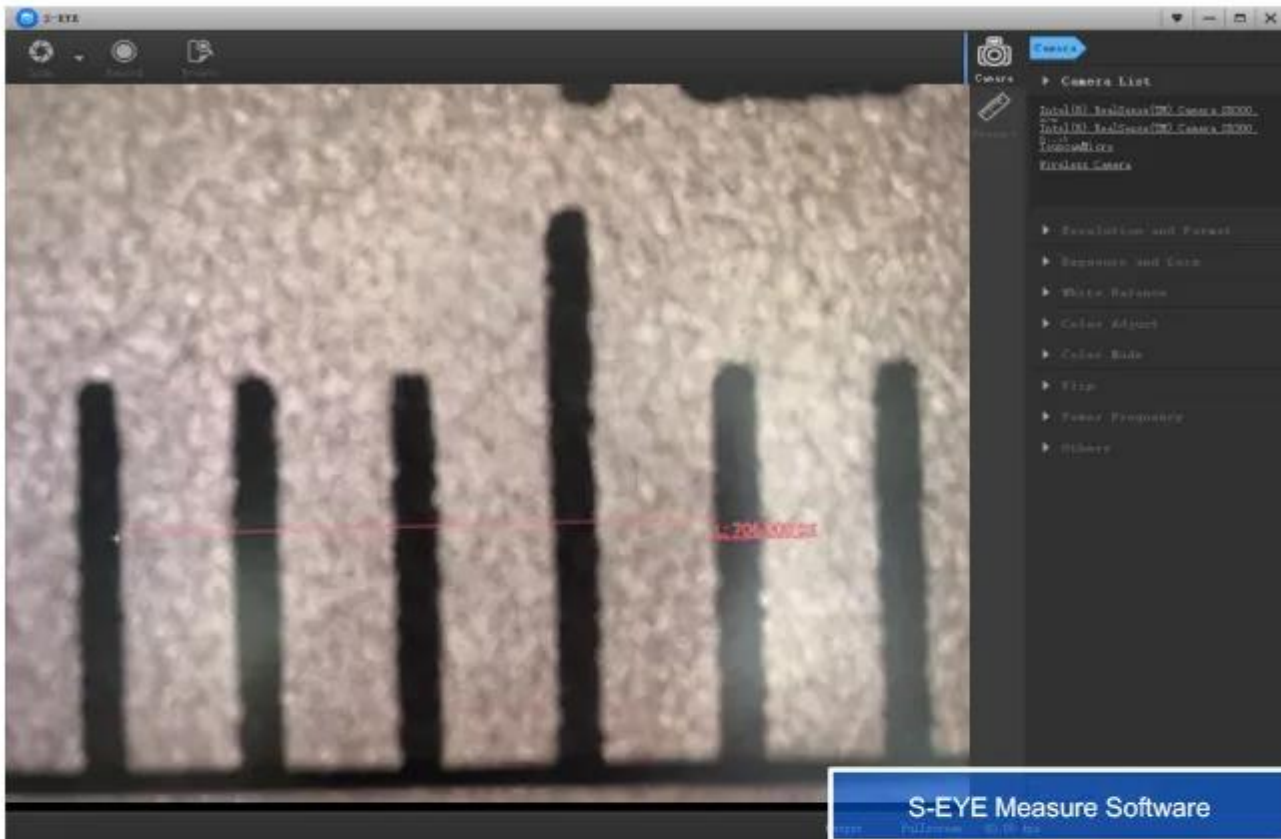

Product Details



TECHNICAL SPECIFICATION

Items	11.6"Win10 Tablet Digital Camera, 16.0M
Model	SLV-16MP TAB
Sensor	1/2.3"Sony CMOS Sensor
Resolution	16M CMOS
LCD Screen	11.6"LCD, Resolution 1920x1080p IPS Screen
CPU	Inter J4105
System	Win10 64Bits
WIFI	2.4G/5G WIFI, 802.11a/b/g/n/ac
Bluetooth	Bluetooth 4.2
Output	USB3.0, HDMI, WIFI, Ethernet, Type-C Power Supply
Pixel Size	1.335x1.335um
Power Supply	DC 12V/2A Power Adapter, Input 100-240V AC
Function	Photo 16M, Video 1080p Full HD
	Adjust: Brightness, Contrast, Saturation, White Balance, Exposure, Flicker Suppression, Sharpness, Gamma, CB/Cr Gain
	UI Interface: Mouse Control, Touch Screen
	Storage Photo, Video, Preview, SD Card Format
Windows Measure	<p>S-EYE Measure Software: Length: Straight Line Length, Broken Line Length, Curve Length, Parallel Line Distance, Point Line Distance, Line Length Geometry:Radius Circle, Two-point Circle, Three-point Circle, Concentric Circle, Radius Circle, Two-point Circle, Three-point Circle Geometric Area:Polygon, Square Remark: Point Coordinates, Crosshairs, Coordinate Systems, Text Notes Image View Professional Measure Software</p>
Size / Weight	Product Size 285x177x12.5mm, Net Weight 0.9kg
Standard Package	LCD Camera, DC12V Power Adapter
Quick Manual	<ol style="list-style-type: none"> 1. First, connect the product to Type-C DC-12V power supply 2. After connecting the power supply, press the power on button, and the blue light of the product will always be on. 3. Double click the S-EYE Measure Software on the desktop of Window 10 4. Click"Open Video Source"in S-EYE software, and choose camera YW1600 6. Measure and adjust as you need.

